

XenoScreen XL YAS

**Accelerated high-sensitivity microplate assay for the
detection of compounds with androgenic activity**

**using *Saccharomyces cerevisiae* strains
with human androgen (hAR) receptor**

Short protocol

**For 96 data points
Art. No. N32-234-Y**

Upon receipt of your XenoScreen XL YAS Assay kit, **make sure that all reagents are stored appropriately (see pg. 4 for storage instructions)**. If components are damaged or if any problems occur, please contact Xenometrix by phone: ++41-61-482-14-34; fax: ++41-61-482-20-72, or e-mail: info@xenometrix.ch

Change log

Date	New version	Changes
23.05.2018	2.0	<ul style="list-style-type: none">• Creation from the XenoScreen XL YES Short Protocol (v. 2.0).

XenoScreen XL YAS Assay Kit

Table of Contents

Number of Data Points	1
Principle of the Test.....	1
Assay Description.....	2
Safety Precautions	2
Warnings	2
Kit Components.....	3
Storage Conditions	4
Required Equipment and Consumables NOT Included with the Kit	5
Yeast Culture Preparation	6
Plate Layout	7
Sample and Controls Preparation.....	8
Preparation of Test Medium	9
Preparation and Dilution of Yeast Cells	10
Transfer of YES and YAS Yeast Cells to Assay Plates.....	10
Incubation of the Assay Plates	11
Addition of Lysis Buffer and Substrate	11
Reading the Assay Plates	12
Calculations and Data Evaluation	13
Validity Criteria	16
Appendix A Storage of Yeast Cells.....	17

Number of Data Points

This test provides a total of 96 data points.

The assay can be used to test 4 samples in duplicates in 8 concentrations in order to obtain dose-response curves for the calculation of androgenic activities, OR to test 32 single-dose samples in duplicates for simplified “yes - no” answers with an estimation of relative activities. This manual describes both versions of the assay which includes all necessary controls, the calculation of EEQ, the determination of the Limit of Detection LoD, and the Limit of Quantification LoQ.

The free Excel calculation workbook available from Xenometrix is based on the plate layouts described in this manual.

Other plate layouts are possible but are not supported by the calculation sheet.

Principle of the Test

The common Baker's or Brewer's yeast (*Saccharomyces cerevisiae*) was genetically modified to identify compounds that can interact with the human androgen receptor hAR. For this purpose the DNA sequences of hAR were stably integrated into the main chromosome of yeast cells. Additionally, the cells also contain an expression plasmid carrying the reporter gene lacZ encoding the enzyme β -galactosidase and androgen (YAS) responsive elements (Routledge, E.J. and Sumpter, J.P. 1996. *Environ. Toxicol. Chem.* **13**:241–8; Sohoni, P. and Sumpter, J.P. 1998. *Endocrinol.* **158**:327–39).

Upon binding of a ligand, the hAR interacts with the response element on the expression plasmid and modulates the transcription of the lacZ reporter gene. The β -galactosidase is released into the medium and converts the yellow substrate chlorophenol red- β -D-galactopyranoside (CPRG) into red product which can be quantified colorimetrically at 570 nm. The measured OD₅₇₀ correlates directly with the amount of secreted β -galactosidase and thus with the activity of the test substance which binds to the androgen receptor.

The XenoScreen XL YAS assay can identify activating (agonistic) activities of test compounds. For the determination of antagonistic activities Xenometrix offers the XenoScreen XL YES/YAS assay (Art. No. N05-233-Y) which detects both androgenic and androgenic activating and inhibiting activities (agonists and antagonists).

The assay can be used for either water samples or for samples dissolved in a solvent like DMSO. Samples dissolved in a solvent have to be diluted 100x in the assay in order to have acceptable levels of solvent.

The XenoScreen XL YAS uses lyticase and a detergent to facilitate the secretion of the intracellularly synthesized β -galactosidase (Schultis T. and Metzger J.W., 2004. *Chemosphere.* **57**:1649–55). This allows to reduce the incubation time from 48 hrs in the standard YAS assay to 18 hours. In addition the accelerated protocol leads also to an enhanced sensitivity.

Assay Description

Growing yeast cells stably transformed with hER α and a β -galactosidase reporter system are exposed to test compound, positive control chemicals (5 α -dihydrotestosterone). The cells are incubated for 18 hrs at 31°C. The induced cells are lysed in the presence of the yellow substrate CPRG which turns purple in the presence of β -galactosidase. Yeast cell growth is assessed prior to addition of the lysis buffer at an OD₆₉₀. The color development is measured at 570 nm and is corrected for diffraction by cells and debris by a simultaneous measurement of OD₆₉₀. The results are evaluated for androgenic activity, as well as for yeast growth inhibition or cytotoxicity.

Safety Precautions

- The YAS yeast cells are genetically modified organisms (GMO). Please consult with your institutional and regulatory authorities for the requirements for handling, storage and disposal of such organisms in accordance with directive 2009/41/EC of the European Parliament and of the Council of 6 May 2009 on the contained use of genetically modified micro-organisms (replaces Council Directive 90/219/EEC of 23 April 1990).
- The control chemical 5 α -dihydrotestosterone provided in this kit is a hormonally active substance. Please consult the Material and Safety Data Sheets (MSDS) for information on handling, disposal and personal protection.
- Not for use in humans and animals. For research purposes only.
- Do not drink, eat, smoke, or apply cosmetics in designated work areas. Wear laboratory coats and gloves when handling specimens and kit reagents. Wash hands thoroughly afterwards. Do not pipette by mouth.

Warnings

» Please observe all highlighted warnings and hints in the text! «

- **Due to the high sensitivity of the XenoScreen XL YAS assay all containers and pipettes coming into contact with the cells or reagents must be absolutely clean and devoid of any residual chemicals such as detergents.**
- **When reusable items are used they should be thoroughly rinsed with distilled water and ethanol (without any additives). We highly recommend to wear gloves also for the handling of glassware and plastic ware.**
- **All solvents should be of the highest available purity grade without any additives.**
- **Read the whole Instructions for Use before starting the assay!**

Kit Components

Each XenoScreen XL YAS assay kit contains media and reagents for the analysis of either 4 test compounds in 8 dilutions in duplicates, or for 32 test compounds at 1 concentration in duplicates for agonistic androgenic endocrine activity. Positive and solvent controls are included.

Use your own aliquots of YAS yeast cells or order them separately: XenoScreen YAS Strain, 1 vial, Art. No. N05-230-A.

Alternative plate layouts, dilution schemes or replicate numbers are possible, but are not described in this manual and are not supported by the Excel calculation sheet provided by Xenometrix.

Kit contents:

- Basal medium
- Vitamin solution
- L-aspartic acid solution
- L-threonine solution
- Cu-sulfate solution
- 10× DO medium
- 10× SD medium
- Streptomycin/Ampicillin solution
- DMSO
- Substrate solution CPRG
- 5α-dihydrotestosterone positive control
- lyticase
- lacZ lysis buffer
- 2-mercaptoethanol
- 3 96-well plates (2 for the assay, 1 for measurement of cell densities)
- 1 gas-permeable plate sealer
- 2 T25 culture flasks with gas-permeable filter caps

To be ordered separately:

- XenoScreen YAS Strain, 1 vial, Art. No. N05-230-A

Storage Conditions

Each Xenometrix XenoScreen XL YAS kit is shipped at ambient temperature. Please contact Xenometrix if you received the kit later than 10 days after the shipment date indicated on the delivery note (phone: ++41-61-482-14-34; fax: ++41-61-482-20-72, or e-mail: info@xenometrix.ch).

The shipment contains the following components which should be stored **immediately upon arrival** as follows:

-20°C:

- Positive control DHT before reconstitution
- Lyticase
- Streptomycin/Ampicillin
- *If ordered separately: Yeast cells on filter discs when stored for more than 1 month*

2 – 8°C

- CPRG substrate solution
- Basal medium
- Vitamin solution
- L-threonine solution
- 2-mercaptoethanol
- 10× DO medium
- 10× SD medium
- Positive control DHT after reconstitution
- *If ordered separately: Yeast cells on filter discs when used within 1 month*

20 – 25°C (room temperature, liquids protected from light):

- L-aspartic acid solution
- Cu(II)-sulfate solution
- DMSO
- lacZ lysis buffer
- 96-well plates
- Culture flasks
- Plate sealer

Required Equipment and Consumables NOT Included with the Kit


- XenoScreen YAS Strain, 1 vial, Art. No. N05-230-E
- Incubator (31°C) with rotating platform (orbital movement)
- High-humidity container (e.g. plastic box with a tight lid, and wet paper tissue)
- Microplate reader capable to read at 570 and 690 nm
- Adjustable micropipettes and sterile tips
(needed volumes: 2, 20–200, 100–1000 µl)
- Adjustable 8-channel pipettes (needed volumes: 20–200 µl)
- Serological pipettes (sterile)
- Pipetting reservoirs (sterile)
- Gloves
- Freezing tubes for the establishment of your own yeast stocks
- Glycerol for freezing cells

Recommended:

- Inverted microscope to inspect yeast cultures


XenoScreen XL – Assay Procedure

Day -4 – Day -2:

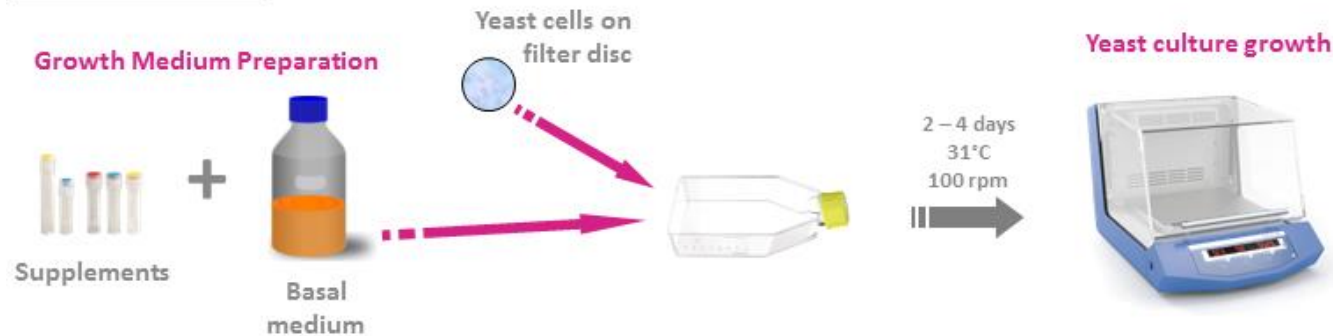


Day 1:

1. Dilution Preparation
directly in assay plate (samples, controls)
2. Addition of yeast cells.

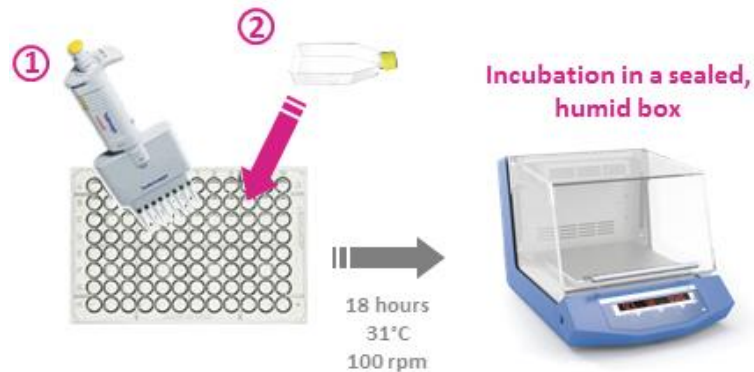

Day 2:

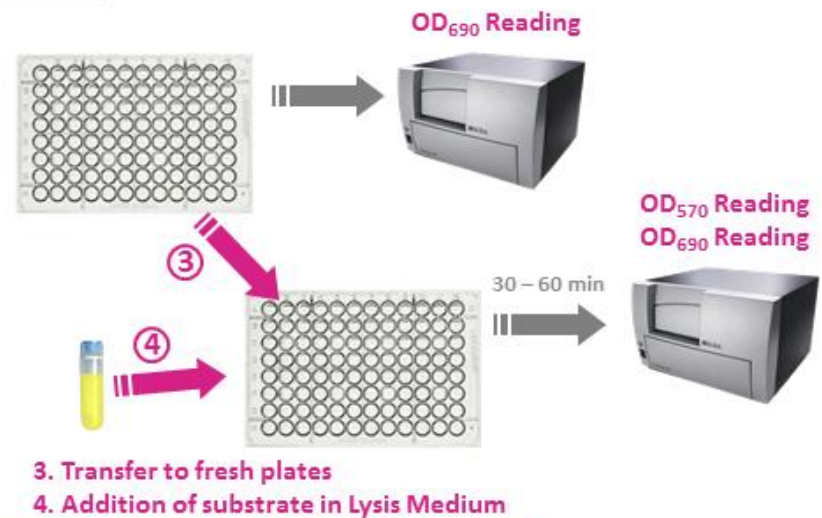

XenoScreen XL

Xenomatrix AG

XenoScreen XL
Assay Procedure