BIOPHEN™ Heparin LRT

REF 221011-RUO R1 R2 4 x 7.5 mL REF 221013-RUO R1 R2 3 x 3 mL REF 221015-RUO R1 R2 4 x 5 mL

Anti-Xa chromogenic method for the assay of Heparin and their analogs, and direct FXa inhibitors, with ready to use liquid reagents.

> FOR RESEARCH USE ONLY. DO NOT USE IN DIAGNOSTIC PROCEDURES.

HYPHEN

155 rue d'Eragny, 95000 Neuville-sur-Oise, France Tél: +33 (0)1 34 40 65 10 Fax: +33 (0)1 34 48 72 36 www.hyphen-biomed.com

English, last revision: 04-2023

info@hyphen-biomed.com

R1 R2 Reagent stability after opening, free from any contamination or evaporation, and stored closed, is of:

- 6 months at 2-8°C.
- 14 days at room temperature (18-25°C).
- Do not freeze.
- Stability on board of the analyzer: see the specific application.

If the substrate becomes yellow, this indicates a contamination. Discard the vial and use a new one.

REAGENTS AND MATERIALS REQUIRED BUT NOT PROVIDED:

- Reagents:

 Distilled water.
- 20% acetic acid or 2% citric acid (endpoint method).
- Physiological Saline (0.9% NaCl).
- Optional: Tris-NaCl-EDTA buffer, pH 7.85 (AR032A/K-RUO), special dilution buffer reducing heparin interferences in direct FXa inhibitors assay
- Specific calibrators and controls with known titration, such as:

	UFH	LMWH	Orgaran [®]	Arixtra [®]
Calibrators	222301-RUO	222001-RUO	222201- RUO	222501- RUO
Controls	223101-RUO / 224101-RUO / 223901-RUO	223001-RUO / 223801-RUO / 224201-RUO / 223701-RUO / 224301-RUO / 224401-RUO	223501- RUO	224001- RUO

	Rivaroxaban standard	Apixaban standard /	Edoxaban standard /	
	/ low range	low range	low range	
0 111 4	222701-RUO /	226201-RUO /	226501-RUO /	
Calibrators	226001-RUO	226101-RUO	226401-RUO	
Controls	224501-RUO /	225301-RUO /	225501-RUO /	
	225101-RUO	225201-RUO	225401-RUO	

Also refer to the specific application guide of the analyzer used

Materials:

- Spectrophotometer or automatic instrument for chromogenic assays.
- Stopwatch; Calibrated pipettes; Silicon glass or plastic test tubes or microplate.

SPECIMEN COLLECTION AND PREPARATION:

The blood (9 volumes) should be carefully collected onto the trisodium citrate anticoagulant (1 volume) (0.109 M, 3.2%) by clean venipuncture. Specific collection

anticoagulant (1 Volume) (0.109 M, 3.2%) by clean Veriptincture. Specific collection tubes for unfractionated heparin testing, such as the CTAD (Citrate, Theophylline, Adenosine and Dipyridamole) tubes, can be used. Discard the first tube. Specimens should be prepared and stored in accordance with applicable local guidelines (for the United States, see the CLSI H21-A5⁶ guideline for further information concerning specimen collection, handling and storage). For plasma storage, please refer to references^{6,7,8,9}.

The kit can be used for kinetics, automated or manual (endpoint) methods. Perform the test at 37°C and read color intensity at 405nm.

For an automated method, application guides are available on request. See specific application guide and specific precautions for each analyzer.

- Assay method (Manual method):

 1. Reconstitute the calibrators and controls as indicated in the specific instructions.
- 2. Dilute the specimens, calibrators (to prepare the calibration curve) and controls in physiological saline buffer, as described in the table below

Products	Calibrators Reference	Controls Reference	Dilution in physiological saline		
LMWH	222001-RUO	223001-RUO / 223801- RUO / 224201-RUO / 223701-RUO / 224301- RUO / 224401-RUO	1:2		
UFH	222301-RUO	223101-RUO / 224101- RUO / 223901-RUO	1:2		
Arixtra®	222501-RUO	224001-RUO	1:2		
Orgaran [®]	222201-RUO	223501-RUO	1:2		
	0.111.4		Dilution in		
Products	Calibrators Reference	Controls Reference	physiological saline or AR032A/K-RUO		
Products Rivaroxaban Standard range		Controls Reference 224501-RUO			
Rivaroxaban Standard range Rivaroxaban Low range	Reference		or AR032A/K-RUO		
Rivaroxaban Standard range	Reference 222701-RUO	224501-RUO	or AR032A/K-RUO 1:10		
Rivaroxaban Standard range Rivaroxaban Low range	222701-RUO 226001-RUO	224501-RUO 225101-RUO	or AR032A/K-RUO 1:10 1:3		
Rivaroxaban Standard range Rivaroxaban Low range Apixaban Standard range	222701-RUO 226001-RUO 226201-RUO	224501-RUO 225101-RUO 225301-RUO	or AR032A/K-RUO 1:10 1:3 1:15		

Establish the calibration curve and test it with the quality controls quickly for optimal assay performance. If stored at room temperature (18-25°C), test the diluted specimens quickly. The exact calibrators and controls concentrations for each batch are indicated on the flyer provided with the kit.

INTENDED USE:

The BIOPHEN™ Heparin LRT kit is an anti-Xa chromogenic method for the *in vitro* quantitative determination of heparin and theirs analogs, in human citrated plasma, quantitative determination of neparin and theirs analogs, in numan citrated piasma, using a manual or automated method. This method is appropriate for the Apixaban, Rivaroxaban and Edoxaban assay, direct Factor Xa (FXa) inhibitors. This method is also appropriate for the determination of anti-Xa activity assay of Arixtra® (Fondaparinux) and Orgaran® (Sodium Danaparoïd), indirect inhibitors whose activity is mediated by plasma antithrombin (AT). Reagents are in the liquid presentation, ready to use (LRT, Liquid reagent Technology).

This kit is for research use only and must not be used for patient diagnosis or treatment.

SUMMARY AND EXPLANATION:

Technical:

Heparin is a sulphated polysaccharide with an high affinity for antithrombin (AT). When complexed with heparin, AT exhibits a fast acting and potent inhibitory activity for coagulant serine esterases: IXa, Xa, XIa, XIIa and thrombin^{1,2}. Low Molecular Weight Heparin (LMWH), and heparin analogues, such as sodium danaparoïd, inhibit more efficiently FXa than thrombin. Anti-Xa assays are then the methods of choice for measuring heparins and their analogues.^{2,5}

The BIOPHEN™ Heparin LRT is a chromogenic anti-Xa method developed for measuring homogeneously unfractionated heparin (UFH) and LMWH, using the same calibration curve. This method is also useful for the determination of anti-Xa activity of Orgarane (sodium danaparoid) and Arixtra (Fondaparinux), indirect inhibitors whose activity is mediated by plasma AT, and for the determination of direct anti-Xa inhibitors (Rivaroxaban, Apixaban and Edoxaban), using specific calibrations.

PRINCIPLE:

The BIOPHEN™ Heparin LRT method is a one stage chromogenic assay based on the inhibition of a constant amount and in excess of FXa, by heparin (or other anti-Xa substance) to be assayed, in the presence of endogenous AT. The residual FXa hydrolyses a specific chromogenic substrate (SXa-11) releasing paranitroaniline (pNA)3. The amount of pNA released (measured by absorbance at 405 nm) is inversely proportional to the concentration of heparin (or other anti-Xa substance) present in the reaction medium.

Heparin and analogues

Heparin + AT → [AT Hep.]

[AT Hep.] + [FXa (excess)] → [FXa-AT-Hep.] + [residual FXa]

[FXa (residual)] + SXa-11 → Peptide + pNA
Rivaroxaban / Apixaban / Edoxaban (DiXal)

[DiXal] + [FXa (excess)] → [FXa-DiXal] + [FXa residual]

[FXa (residual)] + Substrate → Peptide + pNA

R1 Reagent 1: Chromogenic substrate specific for Factor Xa (SXa-11), liquid form. Contains Proclin

R2 Reagent 2: Bovine Factor Xa, liquid form. Contains BSA, Dextran Sulfate4 and small amounts of sodium azide (0.9 g/L).

REF 221011-RUO → R1 R2 4 vials of 7.5 mL.

REF 221013-RUO → R1 R2 3 vials of 3 mL. REF 221015-RUO → R1 R2 4 vials of 5 mL.

WARNINGS AND PRECAUTIONS:

- Some reagents provided in these kits contain materials of animal origin. Users of reagents of these types must exercise extreme care in full compliance with safety precautions in the manipulation of these biological materials as if they were infectious.
- In contact with lead or copper pipes, sodium azide can generate explosive
- Waste should be disposed of in accordance with applicable local regulations.
- Use only the reagents from the same batch of kits.
- Aging studies show that the reagents can be shipped at room temperature without degradation.
- This device of *in vitro* use is intended for professional use in the laboratory

REAGENT PREPARATION:

R1 R2 Reagent is ready to use; homogenize and load it directly on the analyzer following application guide instruction.

For manual method, allow to stabilize for 30 minutes at room temperature (18-25°C), homogenize before use.

STORAGE AND STABILITY:

Unopened reagents should be stored at 2-8°C in their original packaging. Under these conditions, they can be used until the expiry date printed on the kit.

3. Dispense the following to a microplate or plastic tube incubated at 37°C:

	Heparin Orgaran [®] Arixtra [®]		Rivaroxaban Apixaban Edoxaban			
	Microplate	Volume	Volume			
Plasma to test, control or calibrator diluted	30 µL	100 µL	100 μL			
R1 Substrate SXa-11 Preincubated at 37°C	75 µL	250 µL	250 μL			
Mix and incubate at 37°C, for 2 minutes, then introduce:						
R2 Factor Xa Preincubated at 37°C	75 µL	250 µL	250 μL			
Mix and incubate at 37°C for exactly:	60 sec.	60 sec.	120 sec			
Stop the reaction by introducing:						
Citric acid (2%)*	100 μL	350 µL	400 µL			
Mix and measure the absorbance at 405nm against the corresponding blank.						

^{*}Or acetic acid (20%). The yellow color is stable for 2 hours.

The specimen blank is obtained by mixing the reagents in the reverse order to that of the test: Citric acid (2%), R2, R1, diluted specimen.

Measure the optical density at 405 nm. Subtract the measured blank value from the absorbance measured for the corresponding test.

Create a plasma blank if sample is icteric, lipaemic, haemolysed, or if its color differs from the standard plasmas.


The assay of Rivaroxaban, Apixaban and Edoxaban (standard range) can be achieved by kinetics method by recording the change in absorbance between 10 and 35 seconds after the addition of FXa (i.e., Δ OD 405nm). In this case it is not necessary to subtract the blank sample, or to stop the reaction.

If a reaction volume other than that specified above is required for the method used, the ratio of volumes must be strictly observed to guarantee assay performance. The user is responsible for validating any changes and their impact on all results.


CALIBRATION:

The BIOPHEN™ Heparin LRT assay can be calibrated for the assay of various anti-Xa analytes: Heparins and analogs, Rivaroxaban, Apixaban, Edoxaban. Kits containing plasma calibrators specific to these analytes and covering the calibration range are available from HYPHEN BioMed (see the REAGENTS AND MATERIALS REQUIRED BUT NOT PROVIDED paragraph) and can be used to establish the calibration curve specific to the assayed analyte.

<u>Heparin, Arixtra® and Orgaran®:</u> The calibration curves shown below are given by way of example only. The calibration curve established for the assay series must be used.

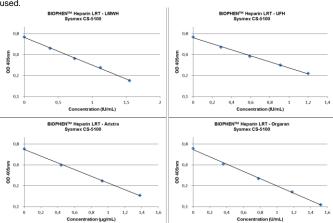

Apixaban, Rivaroxaban, Edoxaban: The calibration curves shown below are given by way of example only. The calibration curve established for the assay series must be used.

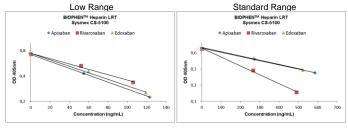

QUALITY CONTROL:

The use of quality controls serves to validate method compliance, along with betweentest assay homogeneity for a given batch of reagents.

Include the quality controls with each series, as per good laboratory practice, in order to validate the test. A new calibration curve should be defined, preferably for each test series, and at least for each new reagent batch, or after analyzer maintenance, or when the measured quality control values fall outside the acceptable range for the method. Each laboratory must define its acceptable ranges and verify the expected performance in its analytical system.

RESULTS:

- For the manual endpoint method, plot the calibration curve, with the OD 405 nm along the Y-axis and the concentration along the X-axis:
 - Rivaroxaban low range, Edoxaban standard range, Arixtra®, Orgaran®, UFH and LMWH, use a Lin-Log scale (concentration-OD).
 - Rivaroxaban standard range, Apixaban and Edoxaban Low range, use a Lin-Lin scale (concentration–OD).

When employing the kinetic method, use \triangle OD 405 instead of OD 405.

- The concentration of heparin (or other anti-Xa molecule) in the test specimen is directly inferred from the calibration curve, if the standard dilution is used.
- If other dilutions are used, the level obtained should be multiplied by the additional dilution factor used.
- Results are expressed in International Units/mL (IU/mL) for Heparin, in U/mL for Orgaran®, in μg/mL for Arixtra®, or in ng/mL for Rivaroxaban, Apixaban and Edoxaban.
- The results obtained should be for research use only and must not be used for patient diagnosis or treatment.

LIMITATIONS:

- To ensure optimum test performance and to meet the specifications, the technical instructions validated by HYPHEN BioMed should be followed carefully.
- Any reagent presenting an unusual appearance or showing signs of contamination must be rejected.
- Any suspicious samples or those showing signs of activation must be rejected.
- Blood activation during collection and plasma preparation, may induce release of Platelet Factor 4 (PF4). PF4 is an inhibitor of heparin. This assay was designed for minimizing the interference of anti-heparin substances in plasma, and especially that of PF4.
- If the AT concentration in the tested plasma is <50%, heparin, Arixtra® or Orgaran®
 can be underestimated as the result of lack of AT (the lack of AT must be confirmed
 by an assay). A variant protocol, with an exogenous source of AT, must then be
 used. High AT concentrations (> 150%) could interfere with the assay.
- High concentration sample may be diluted in a normal plasma pool.
- When a unique curve is used (LMWH/UFH), check that the instrument and application used allow a good superimposition between LMWH and UFH calibrations
- Anti-FXa activity levels may be overestimated in samples containing specific Factor Xa reversal agents, such as andexanet alfa.

PERFORMANCES:

- The lower limit and the measurement range are defined by the analytical system used.
- For the standard range, the calibration range is of about 0 to 600 ng/mL Rivaroxaban/Edoxaban/ Apixaban.
- For the low range, the calibration range is of about 0 to 100 ng/mL Rivaroxaban and of about 0 to 120 ng/mL Apixaban/Edoxaban.
- The calibration range is of about 0 to 1.50 IU/mL for UFH, 0 to 1.75 IU/mL for LMWH, 0 to 1.60 μg/mL for Arixtra® and 0 to 1.75 U/mL for Organan®.
- The enzymatic reaction is rapid, and allows obtaining a high sensitivity for this assay.
- Performance studies were conducted internally using CS-series. Performance was assessed using laboratory controls over a minimum of 10 series and at least one repetition for each control level. The following results were obtained:

	Intra assay			Inter assays				
Sample	n	Mean	CV%	SD	N	Mean	CV%	SD
UFH level 1	10	0.19 IU/mL	2.8	0.01	20	0.20 IU/mL	5.7	0.01
UFH level 2	10	0.56 IU/mL	0.9	0.01	20	0.57 IU/mL	1.3	0.01
LMWH level 3	10	0.78 IU/mL	1.0	0.01	20	0.80 IU/mL	1.1	0.01
LMWH level 4	10	1.22 IU/mL	0.6	0.01	20	1.18 IU/mL	1.2	0.02
Rivaroxaban	30	317 ng/mL	0.9	2.80	20	310 ng/mL	1.2	3.60
Rivaroxaban Low	30	81 ng/mL	0.7	0.54	22	85 ng/mL	4.1	3.51
Apixaban	30	207 ng/mL	1.2	2.53	20	212 ng/mL	2.6	5.58
Apixaban Low	30	85 ng/mL	1.2	0.98	20	84 ng/mL	2.2	1.82
Edoxaban	40	314 ng/mL	0.9	2.82	120	306 ng/mL	1.5	4.68
Edoxaban Low	40	85 ng/mL	1.8	0.70	120	86 ng/mL	3.8	3.22
Orgaran®	10	1.00 U/mL	0.5	0.01	20	1.00 U/mL	0.9	0.01
Arixtra®	10	1.18 µg/mL	0.6	0.01	20	1.19 µg/mL	0.6	0.01

 Correlation with reference method (LCMS:MS vs BIOPHEN™ Heparin LRT, Edoxaban):

Sysmex CS-5100 : n = 144 y = 0.98x + 2.43 r = 0.998

For other molecules, refer to the specific application guide of the analyzer used.

• Interferences: see specific application guide for each analyzer.

REFERENCES:

- Hemker HC and Beguin S. The mode of action of heparin In-vitro and in-vivo. Adv Exp Med Biol (1992).
- Med Biol (1992).

 Leslie B et al. Investigation of the anticoagulant mechanism of a covalent antithrombin-heparin complex J Biol Chem (2000).

 Marlowe CK et al. Design synthesis and structure activity relationship of a series of
- Marlowe CK et al. Design synthesis and structure activity relationship of a series of arginine aldehydes factor Xa Inhibitors. Part 1: structure based on the (D)-Arg-Gly-Arg tripeptide sequence. Bioorg Med Chem Lett (2000).
- Lyon SG et al. Modification of an Amidolytic Heparin Assay to Express Protein-Bound Heparin and to Correct for the Effect of Antithrombin III Concentration. Thromb Hemost (1987).
- Gray E. et al. Heparin and low-molecular-weight heparin. Thromb Haemost. (2008)
 CLSI Document H21-A5: "Collection, transport, and processing of blood specimens for
- CLSI Document H21-As: "Collection, transport, and processing of blood specimens for testing plasma -based coagulation assays and molecular hemostasis assays; approved guideline". (2008).
- Douxfils J, et al. Non-VKA Oral Anticoagulants: Accurate Measurement of Plasma Drug. BioMed Research International. (2015).
- Mauge L. and Alhenc-Gelas M. Stabilité pré-analytique des paramètres de la coagulation: revue des données disponibles. Ann Biol Clin. (2014)
- Birri N et al. Stability of low molecular weight heparin anti-factor Xa activity in citrated whole blood and plasma. Br J Haematol (2011).

SYMBOLS:

Symbols used and signs listed in the ISO 15223-1 standard, see Symbol definitions document.

R1 H317: May cause an allergic skin reaction.

Changes compared to the previous version.