




Manufactured By: HYPHEN BioMed

VITAMIN K DEPENDENT PROTEINS

All the vitamin K dependent factors are extracted from fresh frozen and thawed human citrated plasma with the following procedure: adsorption-elution on baryum citrate; ammonium sulfate salt fractionation; combination of ion exchange, affinity and gel filtration chromatographies.

The final preparation is very homogeneous and keeps the full activity of the native protein. All purification steps are performed with great care and using specific cautions, in order to avoid any activation during the process. No inhibitory activity is present in the final product, which is then stabilized and lyophilized.

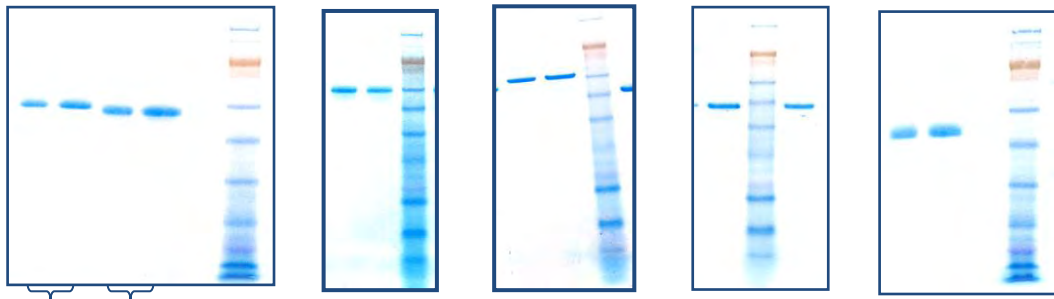

Protein S Factor X

Factor X

Prothrombin

Factor VII

Factor IX



7768 Service Center Drive • West Chester OH 45069

Phone: 513.770.1991

Toll Free: 866.783.3797

Fax: 513.573.9241

Email: info@aniara.com

www.aniara.com

Form AH61
5-2007