

U-CyTech BV

Yalelaan 48  
3584 CM Utrecht  
The Netherlands

P +31.30.253 5960

F +31.30.253 9344

INFO@ucytech.com


www.ucytech.com

## Typical data Marmoset IL-12/23p40 ELISA kit

Catalogue no.: CT345A (5-plate format)

### Validated for:

Marmoset

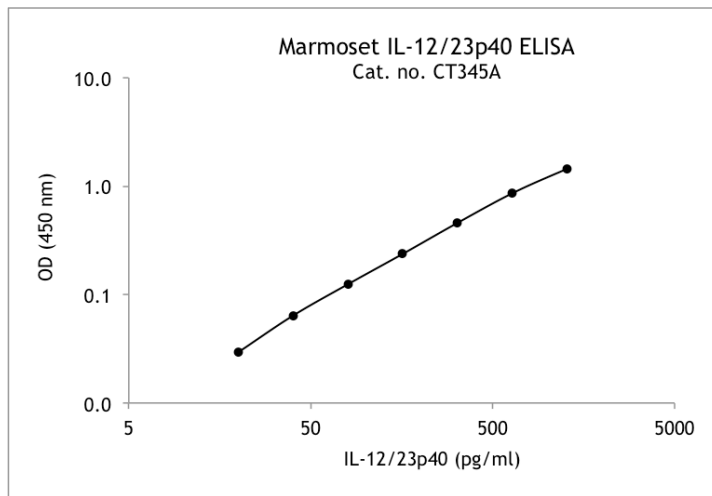

Example of a representative Marmoset IL-12/23p40 standard curve using a 4-parameter logistic fit.

Recommended assay range: 20-1280 pg/ml.

Sensitivity: 20 pg/ml.

Calibration: recombinant and natural marmoset IL-12/23p40.

Type of sample: culture supernatant.

The assay detects free p40 monomers, p40 dimers and the p40 subunits of heterodimeric IL-12 (p35/p40) and heterodimeric IL-23 (p19/p40).

PROUDLY DISTRIBUTED BY

**ANIARA**  
ANIARA DIAGNOSTICA LLC

+1 (513) 770-1991 | +1 (866) 783-3797

7768 Service Center Drive

West Chester, OH 45069

info@aniara.com

Version 140821