


U-CyTech BV
Yalelaan 48
3584 CM Utrecht
The Netherlands
P +31.30.253 5960
F +31.30.253 9344
INFO@ucytech.com
www.ucytech.com

Typical data Human IgG B cell ELISPOT kit


Catalogue no.: CT780-T5 (5-plate format)

Validated for:

Human

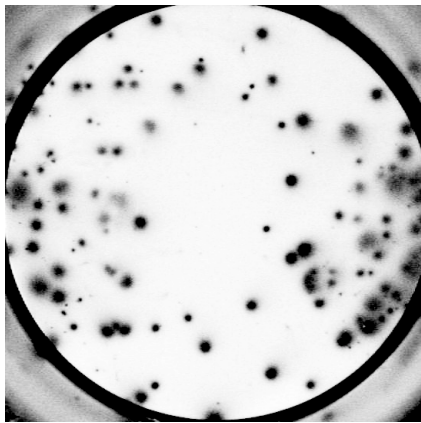

Total number of IgG secreting cells
 2×10^3 cells/well (Assay I)

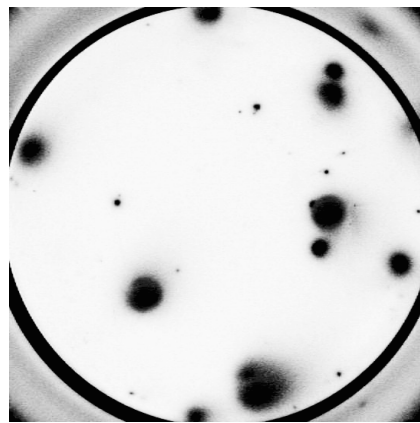

Tetanus toxoid specific IgG secreting cells
 2.5×10^5 cells/well (Assay II)

Example of IgG specific spots produced by human PBMC.

Subject was vaccinated (booster vaccination) intramuscularly with Tetanus toxoid 6 days ago.
B cell activation: 2×10^6 cells/ml, 5 days preincubation in the presence of stimuli.
Incubation ELISPOT: 7 hours.

PROUDLY DISTRIBUTED BY

ANIARA
ANIARA DIAGNOSTICA LLC

+1 (513) 770-1991 | +1 (866) 783-3797
7768 Service Center Drive
West Chester, OH 45069
info@aniara.com

Version 140520