


U-CyTech BV
Yalelaan 48
3584 CM Utrecht
The Netherlands
P +31.30.253 5960
F +31.30.253 9344
INFO@ucytech.com
www.ucytech.com

Typical data Mouse IgG B cell ELISPOT kit


Catalogue no.: CT790-T5 (5-plate format)

Validated for:

Mouse

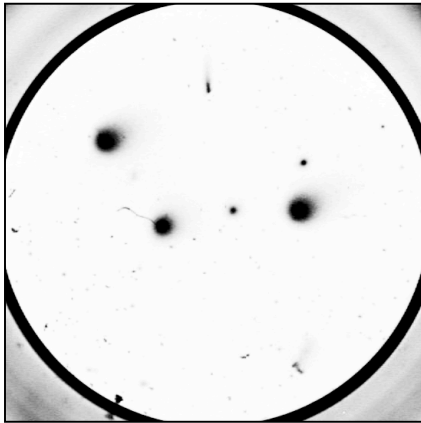

Total number of IgG secreting cells
 1×10^5 cells/well (Assay I)
Stimulus preincubation: none

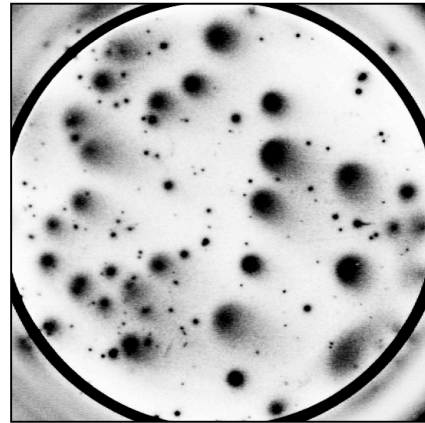

Total number of IgG secreting cells
 2×10^6 cells/well (Assay I)
Stimuli preincubation: IL-2 + R848

Example of IgG specific spots produced by mouse (Balb/c) splenocytes.

B cell activation: 5×10^6 cells/ml, 2 days preincubation without stimulus (left) or in the presence of recombinant IL-2 and R848 (right).

Incubation ELISPOT: 6 hours.

PROUDLY DISTRIBUTED BY

ANIARA
ANIARA DIAGNOSTICA LLC

+1 (513) 770-1991 | +1 (866) 783-3797
7768 Service Center Drive
West Chester, OH 45069
info@aniara.com

Version 140530